



Patient Call Software Driver Installation Guide

www.codasign.co.uk

Once you have chosen which PC you would like the Patient Call box attached to, you will need to install a small piece of software.



Important: You will need admin rights to install the software, so you may need to log a ticket with your IT team for assistance.

1

Download the Driver


Download the driver from the link below and save it to your PC:

<https://downloads.numed.co.uk/Support/EnvisageTools/UC232A.exe>

2

Run the Installer


Open your Downloads folder, Right-click on UC232A_Windows_Setup_v1.0.087 and select 'Run as administrator'.

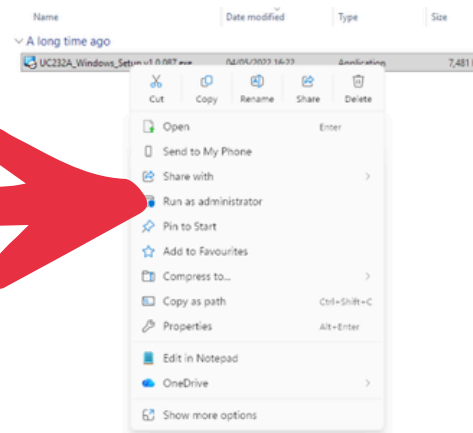

3


Start Installation

Ensure English is selected and click 'Next'.

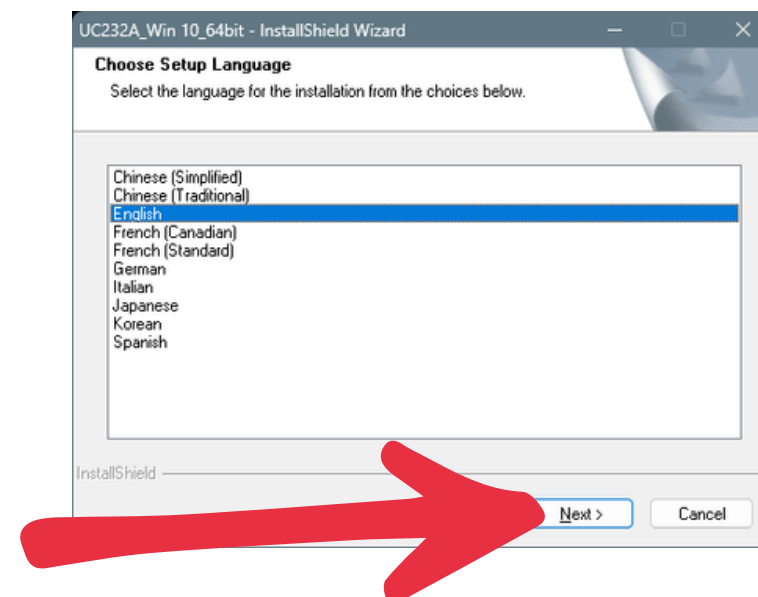


4

Click Next.

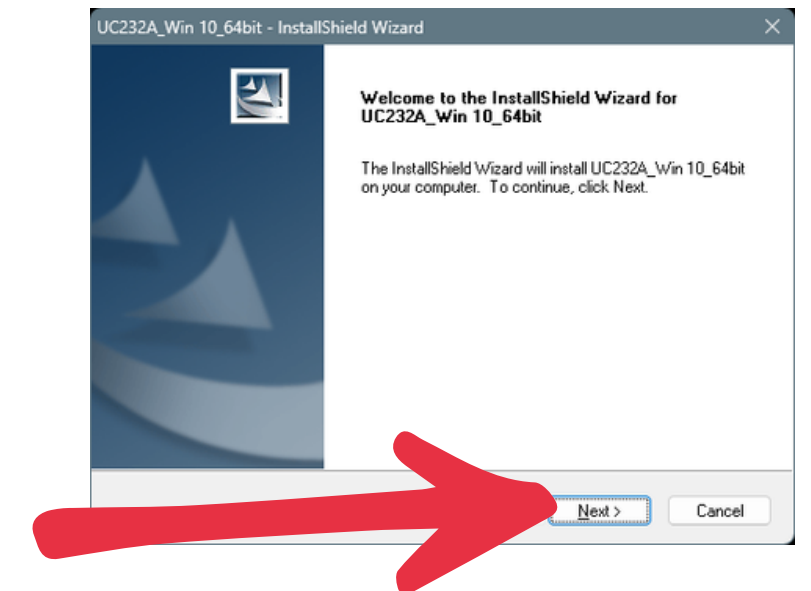

5


Click Install.

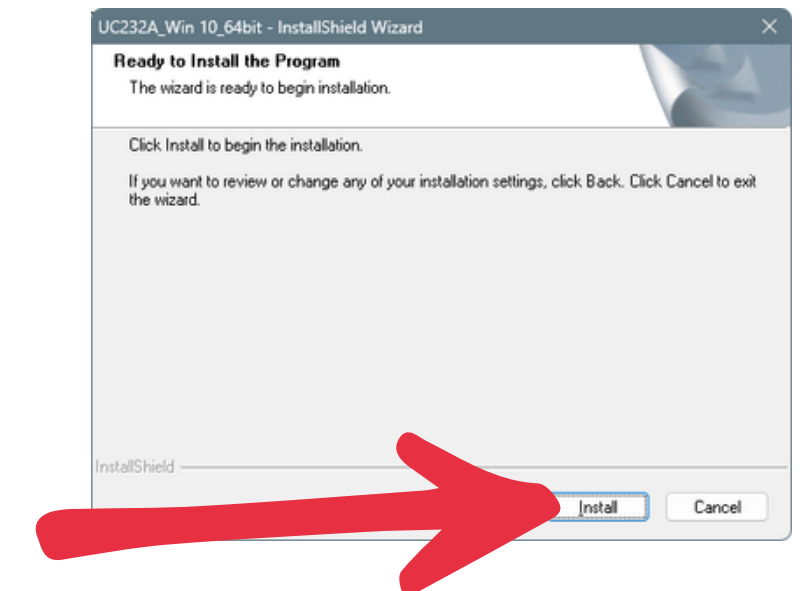

6

Click Finish.

Installation is now complete.

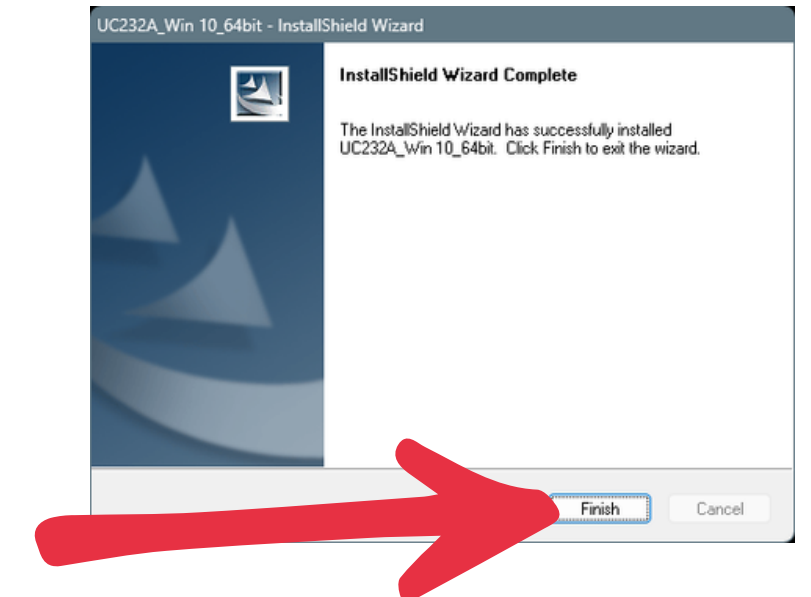



 CodaSign